




Herbicides

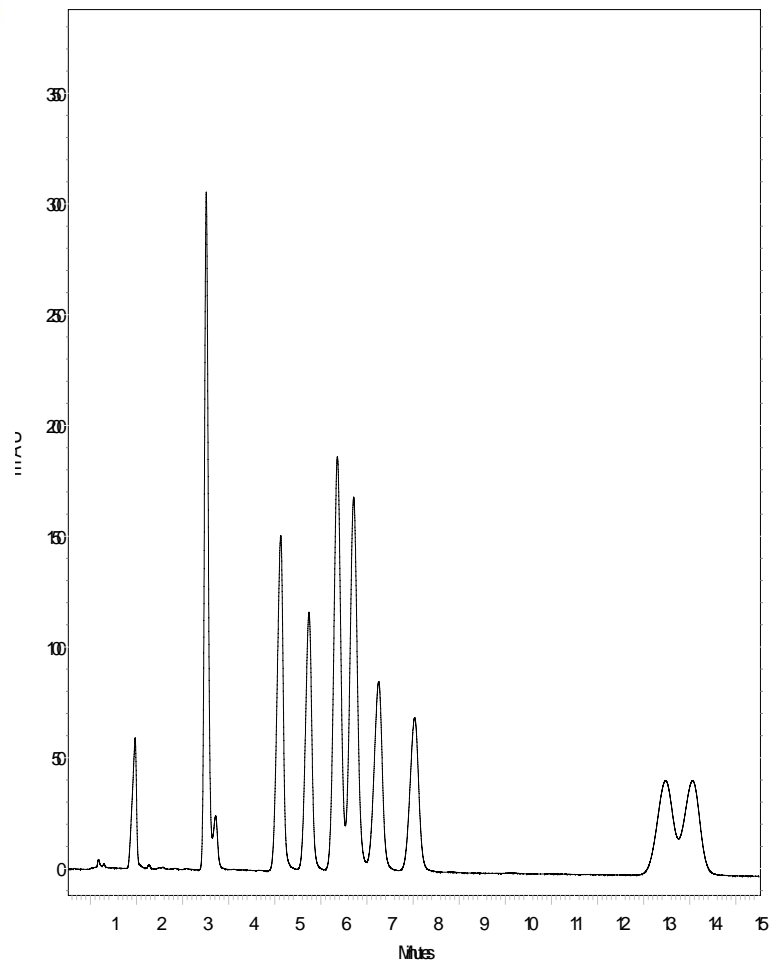

- Column : **Fortis Cyano 50x2.1mm 3 μ**
- p/n : **FCN-020303**
- Mobile Phase:
80:20 – H₂O + 0.2% acetic acid:ACN

- Flow : **0.2ml/min**
- Temp : **25°C**
- λ : **280nm**

- 1. Banzel**
- 2. 2,4-D**
- 3. MCPA**
- 4. PCOC**
- 5. 2,4-DCP**
- 6. 2,4-DP**
- 7. CMPP**
- 8. 2,4-DB**
- 9. MCPB**