




# Positional Isomers

## Hydroxybenzoic acids

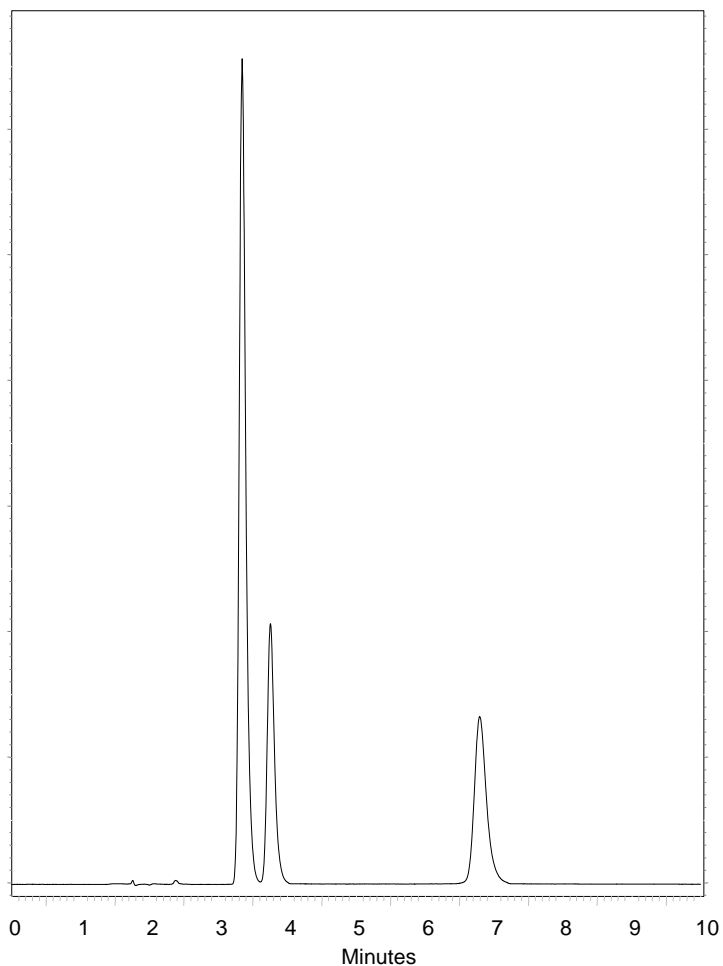

- Column : **Fortis Phenyl 150x4.6mm 5 $\mu$**
- p/n : **FPH-050705**
- Mobile Phase:  
**50:50 – H<sub>2</sub>O + 0.1% formic acid  
:MeOH**
- Flow : **1ml/min**
- Temp : **20°C**
- $\lambda$  : **280nm**

1. **4-Hydroxybenzoic acid**
2. **3-Hydroxybenzoic acid**
3. **2-Hydroxybenzoic acid**

