




Steroids

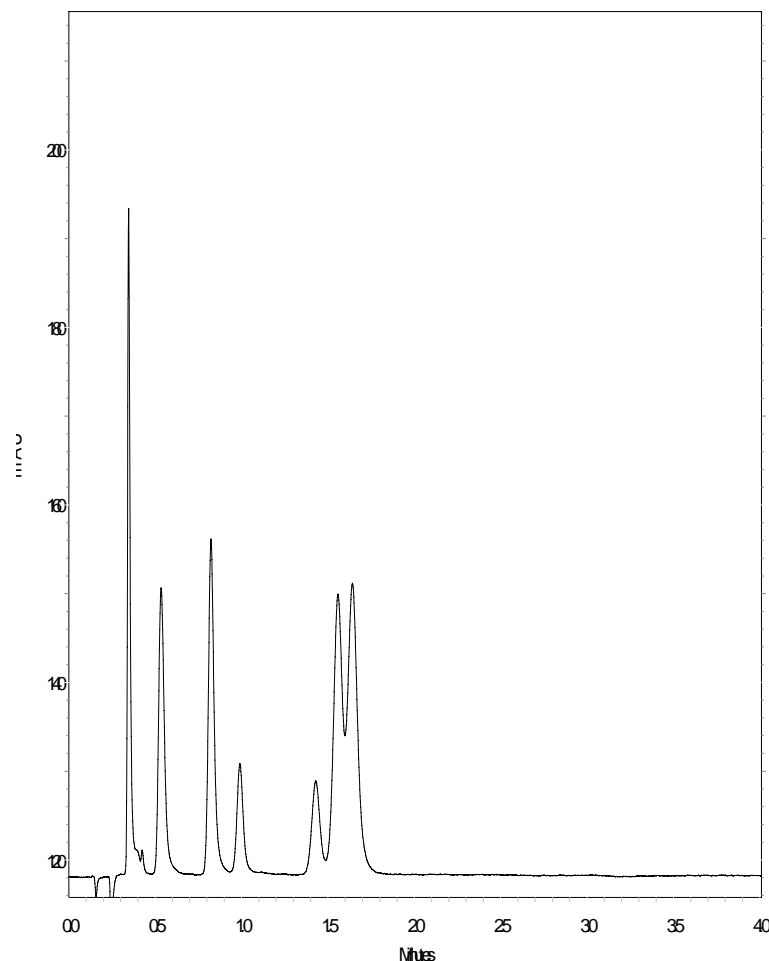

- Column : **Fortis Cyano 50x2.1mm 3 μ**
- p/n : **FCN-020303**
- Mobile Phase:
90:5:5 – Heptane:THF:ACN

- Flow : **0.8ml/min**
- Temp : **25°C**
- λ : **254nm**

- 1. Progesterone**
- 2. 17 α Hydroxyprogesterone**
- 3. 11 α Hydroxyprogesterone**
- 4. Hydroxy-21-acetate**
- 5. Cortisone**
- 6. Prednisone**
- 7. Prednisolone**